

miCON. Because the future belongs to Soundability.

Siemens Life micon at a glance.

- Powered by micon – the new platform behind BestSound Technology
- e2e wireless 2.0
- Battery door on/off
- Programmable push button
- Telecoil and audio input via miniTek™
- Color exchange concept
- Nanocoated housing
- Size 312 battery
- Range (2ccm Coupler): LifeTube 47/121 dB Earhook 57/123 dB



Available in three performance levels

	7mi	5mi	3mi
High res. adaptive TwinMic™ System	48	32	24
Directional speech enhancement	●	–	–
Situation equalizer/no. of situations	6	3	1
High res. SpeechFocus	●	●	–
Frequency compression	●	●	●

Stylish looks with everything inside.

Nanocoating
Unique protection for BTE models from moisture and debris.

Exchangeable housing
The housing can easily be exchanged if the customer's color preferences change.

Push button
Programmable push button acts as a volume control or also allows you to select different hearing programs.

LifeTube
Ultra thin and transparent, the inconspicuous tubing delicately hugs the ear for the most comfortable fit.

Tinnitus noiser
Option for easy, controllable and effective relief.

Battery compartment
With on/off function.


Semi-open tip
For feedback-free semi-open fittings appropriate for a wide range of ear canals.

Chic designs to keep up with the times.

Color Conversion Kit

Siemens Life micon follows all trends. Whenever your customer changes their color preference, you can exchange the complete housing with the help of the Color Conversion Kit. In the kit, you will find all the tools necessary for an easy exchange.

Straightforward housing exchange

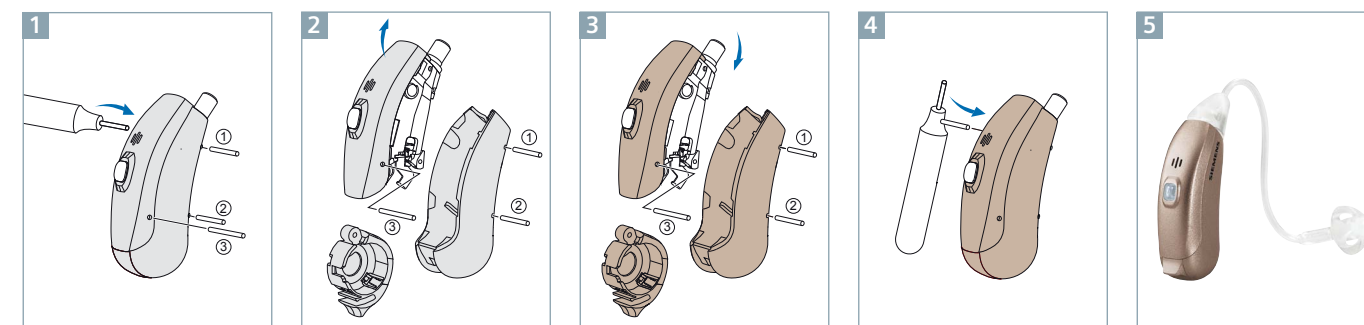


1. Undo the pins using the removal tool.
2. Detach the old housing.
3. Put on the new housing.
4. Replace the pins using the removal tool.
5. And the job is done!


Accessories

ePen, ProPocket and easyPocket

Remote controls for convenient and discreet adjustments. The new easyPocket is extremely ergonomic, has an easy-to-read display and modern, stylish design.

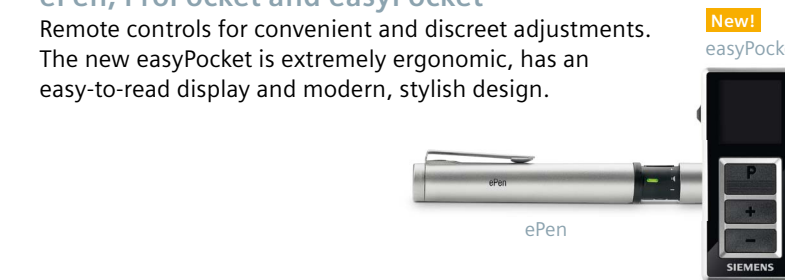

miniTek Remote App

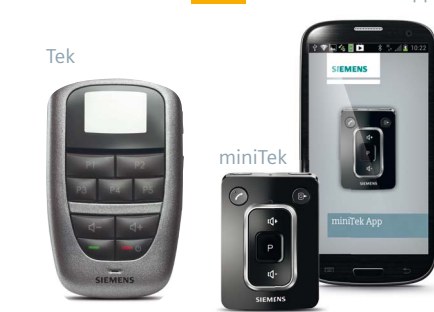

Tek, miniTek and miniTek Remote App

Siemens Tek™ and miniTek work seamlessly with Siemens Life micon, wirelessly connecting the hearing instrument to many different devices. Sound from Bluetooth phones, MP3 players, television and other audio sources is streamed directly to the hearing instruments. And now Android smartphones can be used conveniently to operate the miniTek via the new miniTek Remote App, available in the Google Play* Store.

*Google Play is a trademark of Google Inc.

VoiceLink

The VoiceLink™ companion microphone connects to miniTek and is a practical aid when wearers are out to dinner in a noisy restaurant or attending conferences. VoiceLink can be given to the speaker for direct speech transmission to the user's instruments.




Highlights

- Enhanced tinnitus noiser: option for easy, controllable and effective relief
- Programmable push button, also as volume control
- Semi-open tip: for feedback-free semi-open fittings, appropriate for a wide range of ear canals



Semi-open tips

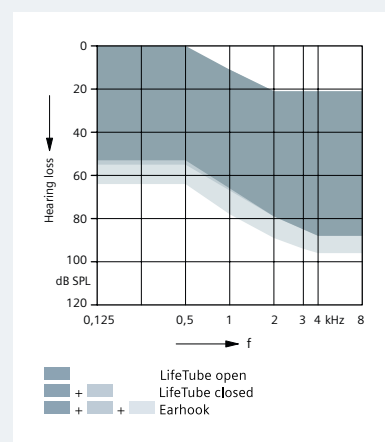

Personalized tinnitus control


Tinnitus noiser

One of the key factors for tinnitus treatment is learning to control the tinnitus – rather than being controlled by it. Many tinnitus treatments such as habituation or masking involve the use of an external acoustic training device. The improved tinnitus noiser of Life micon is an effective control instrument that distracts users from their tinnitus. Not only can the shape of the noise signal be configured in up to 20 bands – this customized masking signal can also be mixed with the microphone signal, providing customized relief in any situation.

The tinnitus function in detail:

- Separate noiser signal generator
- Four pre-programmed therapy signals: white noise, pink noise, speech noise and high tone noise
- Individual fine-tuning of noise program for up to 20 bands
- Three operating modes: microphone signal only, tinnitus noiser function only, mixed mode

Siemens Life micon colors

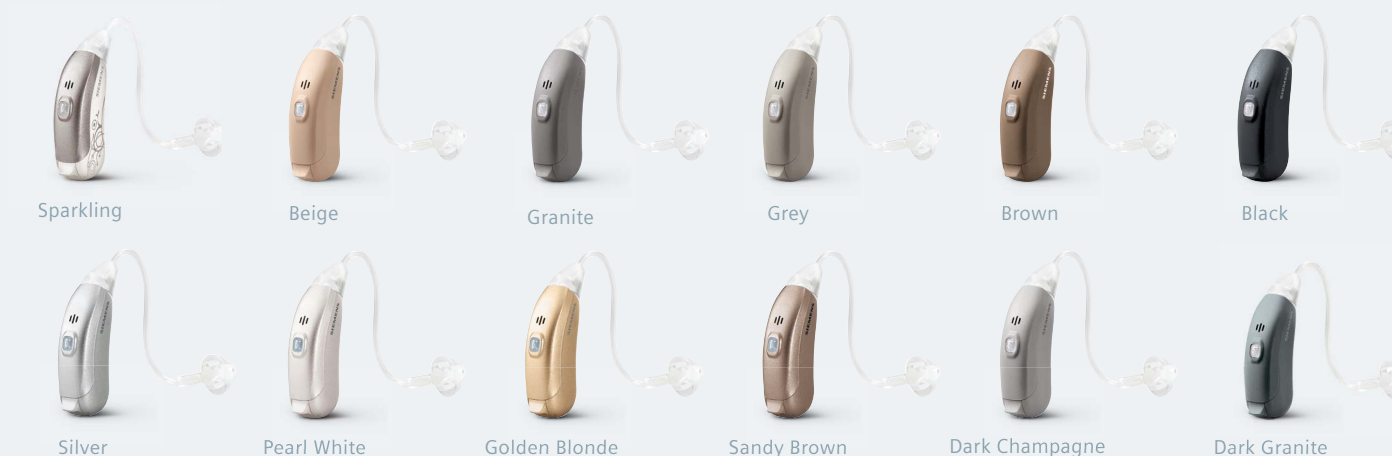

Fitting accessories

Aside from a full range of proven fitting accessories, the new semi-open tips are now also available for feedback-free semi-open fittings.



LifeTips

Open		4, 6, 8, 10 mm
Semi-open		8–12 mm
Closed		8 mm
Double		8/10, 10/12 mm



miCON. Because the future belongs to Soundability.

micon is the new platform behind BestSound Technology. Its three key features combine to provide the perfect balance between Sound Quality and effective Audibility. We call this state of equilibrium Soundability™. The result: maximum spontaneous acceptance, long-term listening satisfaction and a natural sound experience that fits the individual's needs.

miSound

This technology is so transparent that the wearer can forget that they are hearing through a hearing instrument. New micon fit fitting formula, enhanced feedback cancellation and bandwidth extension offer the most natural comfortable sound experience.

miFocus

This new 48-channel TwinMic System™ makes hearing effortless again. Directional microphones, improved localized noise suppression, higher frequency resolution and optional frequency compression for appropriate hearing losses offer more speech intelligibility, even in the most difficult environments.

miGuide

This powerful engine offers the wearer customized sound with automatic situation detection, an intelligent algorithm that learns user sound preferences, and automatic gain adjustment, which increases the gain over time to compensate for the effects of auditory deprivation. With Connexx™ 7, frequency shaping for the detected situations can be manually adjusted.

Creating the balance: BestSound Technology



Soundability – the ideal equilibrium for personalized perfection.

miCON
Because the future belongs to Soundability.



Siemens Life micon

The Bluetooth word mark and logos are owned by the Bluetooth SIG, Inc. and any use of such marks by Siemens Audiologische Technik is under license. Other trademarks and trade names are those of their respective owners.

The information in this document contains general descriptions of the technical options available, which do not always have to be present in individual cases and are subject to change without prior notice.

Global Siemens Headquarters
Siemens AG
Wittelsbacherplatz 2
80333 Muenchen
Germany

Global Siemens Healthcare Headquarters
Siemens AG
Healthcare Sector
Henkestrasse 127
91052 Erlangen
Germany
Phone: +49 9131 84-0
www.siemens.com/healthcare

Legal Manufacturer
Siemens Audiologische Technik GmbH
Gebbertstrasse 125
DE-91058 Erlangen
Germany

Local Contact Information
Siemens contact
House number and street name
City and DE-postal code
Country
Phone: +XX-XXX-XXX-XXXX
www.XXXXXXXXXXXXXXXXXX

Global Business Unit
Siemens Audiologische Technik GmbH
Gebbertstrasse 125
DE-91058 Erlangen
Germany
Phone: +49 9131 308-0

Order No. A91SAT-01829-99C2-7600 | Printed in Germany | © 03.2013 Siemens AG

www.siemens.com/hearing

SIEMENS



www.siemens.com/life-hearing

Siemens Life micon.

Information for Hearing Care Professionals.



Life sounds brilliant.

Sense of hearing meets sense of style.

The elegant Siemens Life features micon™, the new platform behind BestSound™ Technology, to offer the best possible sound comfort from the very start. And it now comes in a new performance level of 3mi. As the smallest standard BTE in the Siemens product range, it is the ideal alternative if a RIC hearing instrument is not an option. Siemens Life micon is appropriate for a wide variety of ear canals, whether bent or simply extremely small, and also if a patient has cerumen issues. Fitting is simple thanks to the flexible LifeTubes – Siemens Life micon fits virtually every user right away. It is simple, convenient and reliable – also thanks to its proven housing.



Siemens Life micon



Siemens researchers won the "Deutscher Zukunftspreis 2012" (German Future Award) for the development of its binaural hearing system. Siemens Life micon offers this technical innovation for improved spatial perception with its e2e wireless™ 2.0 technology.